

# AssayMax™ Human Annexin A5 ELISA Kit

Assaypro LLC 3400 Harry S Truman Blvd St. Charles, MO 63301 T (636) 447-9175 F (636) 395-7419 www.assaypro.com

For any questions regarding troubleshooting or performing the assay, please contact our support team at support@assaypro.com.

Thank you for choosing Assaypro.

## **Assay Summary**

**Step 1**. Add 50  $\mu$ l of Standard or Sample per well. Incubate 2 hours.

Step 2. Wash, then add 50  $\mu l$  of Biotinylated Antibody per well. Incubate 1 hour.

Step 3. Wash, then add 50  $\mu l$  of SP Conjugate per well. Incubate 30 minutes.

**Step 4.** Wash, then add 50  $\mu$ l of Chromogen Substrate per well. Incubate 15 minutes.

**Step 5.** Add 50  $\mu$ l of Stop Solution per well. Read at 450 nm immediately.

# Symbol Key

Consult instructions for use.

# Assay Template

12								
11								
10								
6								
8								
7								
9								
ß								
4								
m								
2								
1								
	A	B	С	D	Е	ł	Ð	н

# AssayMax<sup>™</sup> Human Annexin A5 (ANXA5) ELISA Kit

Catalog No. EA3601-1 Sample insert for reference use only

#### Introduction

Annexin A5 (ANXA5, Annexin V) is a calcium dependent phospholipid binding protein of the annexin group. Although little is known about the function of this protein, it has been proposed that calcium-specific adsorption of annexin A5 to phospholipid monolayers induces polymerization of the protein in the plane of adsorption, potentially exerting a biological effect on the membrane itself or upon membrane-bound proteins. Electron image analysis has shown that annexin A5 molecules form trimers in a triskelion pattern on lipid monolayers (1). Annexin A5 also displays upregulation and has been shown to interact with kinase insert domain receptor and integrin beta 5 (2-4).

#### Principle of the Assay

The AssayMax<sup>™</sup> Human Annexin A5 ELISA (Enzyme-Linked Immunosorbent Assay) Kit is designed for detection of annexin A5 in human **plasma, serum, saliva, cell lysate, and tissue samples**. This assay employs a quantitative **sandwich enzyme immunoassay** technique that measures human annexin A5 in approximately 4 hours. A polyclonal antibody specific for human annexin A5 has been pre-coated onto a 96-well microplate with removable strips. Annexin A5 in standards and samples is sandwiched by the immobilized antibody and a biotinylated polyclonal antibody specific for human annexin A5, which is recognized by a streptavidin-peroxidase (SP) conjugate. All unbound material is washed away and a peroxidase enzyme substrate is added. The color development is stopped and the intensity of the color is measured.

#### **Caution and Warning**

- This product is for **Research Use Only** and is not intended for use in diagnostic procedures.
- Prepare all reagents (diluent buffer, wash buffer, standard, biotinylated antibody, and SP conjugate), as instructed, prior to running the assay.
- Prepare all samples prior to running the assay. The dilution factors for the samples are suggested in this insert. However, the user should determine the optimal dilution factor.

- Spin down the SP conjugate vial, the biotinylated antibody vial, and the standard diluent vial before opening and using contents.
- The Stop Solution is an acidic solution.
- The kit should not be used beyond the expiration date.

#### Reagents

- Human Annexin A5 Microplate: A 96-well polystyrene microplate (12 strips of 8 wells) coated with a polyclonal antibody against human annexin A5.
- **Sealing Tapes:** Each kit contains 3 precut, pressure sensitive sealing tapes that can be cut to fit the format of the individual assay.
- Human Annexin A5 Standard: Human annexin A5 in a buffered protein base (192 ng, lyophilized).
- **Biotinylated Human Annexin A5 Antibody (50x):** A 50-fold concentrated biotinylated polyclonal antibody against human annexin A5 (120 µl).
- **EIA Diluent Concentrate (10x):** A 10-fold concentrated buffered protein base (20 ml).
- Standard Diluent (1x): A buffered protein base with stabilizer (2 ml).
- Wash Buffer Concentrate (20x): A 20-fold concentrated buffered surfactant (30 ml, 2 bottles).
- SP Conjugate (100x): A 100-fold concentrate (80 µl).
- **Chromogen Substrate (1x):** A stabilized peroxidase chromogen substrate tetramethylbenzidine (7 ml).
- Stop Solution (1x): A 0.5 N hydrochloric acid solution to stop the chromogen substrate reaction (11 ml).

## **Storage Condition**

- Upon arrival, immediately store components of the kit at recommended temperatures up to the expiration date.
- Store Standard, SP Conjugate, and Biotinylated Antibody at -20°C.
- Store Microplate, Diluent Concentrate (10x), Standard Diluent (1x), Wash Buffer, Stop Solution, and Chromogen Substrate at 2-8°C.
- Unused microplate wells may be returned to the foil pouch with the desiccant packs and resealed. May be stored for up to 30 days in a vacuum desiccator.

## **Other Supplies Required**

- Microplate reader capable of measuring absorbance at 450 nm
- Pipettes (1-20 μl, 20-200 μl, 200-1000 μl, and multiple channel)
- Deionized or distilled reagent grade water

#### Sample Collection, Preparation, and Storage

- **Plasma:** Collect plasma using one-tenth volume of 0.1 M sodium citrate as an anticoagulant. Centrifuge samples at 3000 x g for 10 minutes and collect plasma. The sample is suggested for use at 1x; however, user should determine optimal dilution factor depending on application needs. The undiluted samples can be stored at -20°C or below for up to 3 months. Avoid repeated freeze-thaw cycles (EDTA or Heparin can also be used as an anticoagulant).
- Serum: Samples should be collected into a serum separator tube. After clot formation, centrifuge samples at 3000 x g for 10 minutes and remove serum. The sample is suggested for use at 1x; however, user should determine optimal dilution factor depending on application needs. The undiluted samples can be stored at -20°C or below for up to 3 months. Avoid repeated freeze-thaw cycles.
- Saliva: Collect saliva using sample tube. Centrifuge samples at 800 x g for 10 minutes. The sample is suggested for use at 1x or within the range of 2x 10x into EIA Diluent; however, user should determine optimal dilution factor depending on application needs. The undiluted samples can be stored at -20°C or below for up to 3 months. Avoid repeated freeze-thaw cycles.
- Cell Lysate: Rinse cell with cold PBS and then scrape the cell into a tube with 5 ml of cold PBS and 0.5 M EDTA. Centrifuge suspension at 1500 rpm for 10 minutes at 4°C and aspirate supernatant. Resuspend pellet in ice-cold Lysis Buffer (PBS, 1% Triton X-100, protease inhibitor cocktail). For every 1 x 10<sup>6</sup> cells, add approximately 100 µl of ice-cold Lysis Buffer. Incubate on ice for 60 minutes. Centrifuge at 13000 rpm for 30 minutes at 4°C and collect supernatant. If necessary, dilute samples into EIA Diluent; user should determine optimal dilution factor depending on application needs. The undiluted samples can be stored at -80°C. Avoid repeated freeze-thaw cycles.
- **Tissue:** Extract tissue samples with 0.1 M phosphate-buffered saline (pH 7.4) containing 1% Triton X-100 and centrifuge at 14000 x g for 20 minutes. Collect the supernatant and measure the protein concentration. If necessary, dilute samples into EIA Diluent; user should determine optimal dilution factor depending on application needs. Store remaining extract at -80°C. Avoid repeated freeze-thaw cycles.

Applicable samples may also include biofluids, cell culture, and tissue homogenates. If necessary, user should determine optimal dilution factor depending on application needs.

#### Refer to Dilution Guidelines for further instruction.

	<b>Guidelines for Dilutions of 100-fold or Greater</b> (for reference only; please follow the insert for specific dilution suggested)			
	100x	10000x		
-				
A)	4 μl sample : 396 μl buffer (100x)	A)	4 μl sample : 396 μl buffer (100x)	
	= 100-fold dilution	B)	4 μl of A : 396 μl buffer (100x)	
			= 10000-fold dilution	
	Assuming the needed volume is less than		Assuming the needed volume is less than	
	or equal to 400 μl.		or equal to 400 μl.	
	1000x		100000x	
A)	4 μl sample : 396 μl buffer (100x)	A)	4 μl sample : 396 μl buffer (100x)	
B)	24 μl of A : 216 μl buffer (10x)	B)	4 μl of A : 396 μl buffer (100x)	
	= 1000-fold dilution	C)	24 μl of B : 216 μl buffer (10x)	
			= 100000-fold dilution	
	Assuming the needed volume is less than		Assuming the needed volume is less than	
	or equal to 240 μl.		or equal to 240 μl.	

#### **Reagent Preparation**

- Freshly dilute all reagents and bring all reagents to room temperature before use.
- **EIA Diluent Concentrate (10x):** Dilute the EIA Diluent Concentrate 10fold with reagent grade water to produce a 1x solution. When diluting the concentrate, make sure to rinse the bottle thoroughly to extract any precipitates left in the bottle. Mix the 1x solution gently until the crystals have completely dissolved. Store for up to 30 days at 2-8°C.
- Human Annexin A5 Standard: Reconstitute the Human Annexin A5 Standard (192 ng) with 0.75 ml of Standard Diluent to generate a 256 ng/ml standard stock solution. Allow the vial to sit for 10 minutes with gentle agitation prior to making dilutions. From the standard stock solution (256 ng/ml), dilute 8-fold with EIA Diluent to produce a 32 ng/ml standard working solution. Prepare duplicate or triplicate standard points by serially diluting the standard working solution (32 ng/ml) 2-fold with equal volume of EIA Diluent to produce 16, 8, 4, 2, 1, and 0.5 ng/ml solutions. EIA Diluent serves as the zero standard (0 ng/ml). Aliquot remaining stock solution to limit repeated freeze-thaw cycles. This solution should be stored at -20°C and used within 30 days.

Standard Point	Dilution	[ANXA5] (ng/ml)
P1	1 part Standard (256 ng/ml) + 7 parts EIA Diluent	32
P2	1 part P1 + 1 part EIA Diluent	16
Р3	1 part P2 + 1 part EIA Diluent	8.0
P4	1 part P3 + 1 part EIA Diluent	4.0
P5	1 part P4 + 1 part EIA Diluent	2.0
P6	1 part P5 + 1 part EIA Diluent	1.0
P7	1 part P6 + 1 part EIA Diluent	0.5
P8	EIA Diluent	0.0

- Biotinylated Human Annexin A5 Antibody (50x): Spin down the antibody briefly and dilute the desired amount of the antibody 50-fold with EIA Diluent to produce a 1x solution. The undiluted antibody should be stored at -20°C.
- Wash Buffer Concentrate (20x): Dilute the Wash Buffer Concentrate 20fold with reagent grade water to produce a 1x solution. When diluting the concentrate, make sure to rinse the bottle thoroughly to extract any precipitates left in the bottle. Mix the 1x solution gently until the crystals have completely dissolved.
- SP Conjugate (100x): Spin down the SP Conjugate briefly and dilute the desired amount of the conjugate 100-fold with EIA Diluent to produce a 1x solution. The undiluted conjugate should be stored at -20°C.

#### **Assay Procedure**

- Prepare all reagents, standard solutions, and samples as instructed. Bring all reagents to room temperature before use. The assay is performed at room temperature (20-25°C).
- Remove excess microplate strips from the plate frame and return them immediately to the foil pouch with desiccants inside. Reseal the pouch securely to minimize exposure to water vapor and store in a vacuum desiccator.
- Add 50 μl of Human Annexin A5 Standard or sample to each well. Gently tap plate to thoroughly coat the wells. Break any bubbles that may have formed. Cover wells with a sealing tape and incubate for 2 hours. Start the timer after the last addition.
- Wash the microplate manually or automatically using a microplate washer. Invert the plate and decant the contents; hit 4-5 times on absorbent material to completely remove the liquid. If washing manually, wash five times with 200 µl of Wash Buffer per well. Invert the plate each time and decant the contents; hit 4-5 times on absorbent material to completely remove the liquid. If using a microplate washer,

wash six times with 300  $\mu$ l of Wash Buffer per well; invert the plate and hit 4-5 times on absorbent material to completely remove the liquid.

- Add 50 µl of Biotinylated Human Annexin A5 Antibody to each well. Gently tap plate to thoroughly coat the wells. Break any bubbles that may have formed. Cover wells with a sealing tape and incubate for 1 hour.
- Wash the microplate as described above.
- Add 50 µl of SP Conjugate to each well. Gently tap plate to thoroughly coat the wells. Break any bubbles that may have formed. Cover wells with a sealing tape and incubate for 30 minutes. Turn on the microplate reader and set up the program in advance.
- Wash the microplate as described above.
- Add 50 μl of Chromogen Substrate to each well. Gently tap plate to thoroughly coat the wells. Break any bubbles that may have formed. Incubate in ambient light for 15 minutes or until the optimal blue color density develops.
- Add 50 µl of Stop Solution to each well. The color will change from blue to yellow. Gently tap plate to ensure thorough mixing. Break any bubbles that may have formed.
- Read the absorbance on a microplate reader at a wavelength of 450 nm **immediately**. If wavelength correction is available, subtract readings at 570 nm from those at 450 nm to correct optical imperfections. Otherwise, read the plate at 450 nm only. Please note that some unstable black particles may be generated at high concentration points after stopping the reaction for about 10 minutes, which will reduce the readings.

#### Data Analysis

- Calculate the mean value of the duplicate or triplicate readings for each standard and sample.
- To generate a standard curve, plot the graph using the standard concentrations on the x-axis and the corresponding mean 450 nm absorbance (OD) on the y-axis. The best fit line can be determined by regression analysis using log-log or four-parameter logistic curve fit.
- Determine the unknown sample concentration from the Standard Curve and multiply the value by the dilution factor.


#### **Typical Data**

• The typical data is provided for reference only. Individual laboratory means may vary from the values listed. Variations between laboratories may be caused by technique differences.

Standard Point	ng/ml	OD	Average OD
P1	32	1.889	1.934
		1.979	
P2	16	1.665	1.672
		1.679	1.072
Р3	8.0	1.219	1.243
FJ	8.0	1.267	1.245
P4	4.0	0.785	0.808
г4	4.0	0.831	0.808
P5	2.0	0.467	0.482
22		0.497	0.482
P6	1.0	0.302	0.205
PO		0.308	0.305
Р7	0.5	0.203	0.198
۲/	0.5	0.193	0.196
DQ.	0.0	0.082	0.084
P8	0.0	0.086	0.084
Sample: Poo	oled Normal	0.428	0.400
Sodium Citrat		0.390	0.409
Sample: Poo	oled Normal	0.623	0.500
Serun	n (1x)	0.569	0.596

#### **Standard Curve**

• The curve is provided for illustration only. A standard curve should be generated each time the assay is performed.

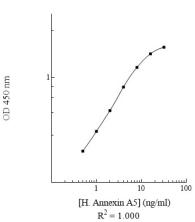

Human Annexin A5 Standard Curve

#### **Reference Value**

• Plasma and serum samples from healthy adults were tested (n=40). On average, human annexin A5 level was 2.06 ng/ml.

Sample	n	Average Value (ng/ml)
Pooled Normal Plasma	10	1.58
Normal Plasma	20	1.64
Pooled Normal Serum	10	2.97

#### **Performance Characteristics**

- This assay recognizes both natural and recombinant human annexin A5.
- The minimum detectable dose of human annexin A5 as calculated by 2SD from the mean of a zero standard was established to be 0.22 ng/ml.
- Intra-assay precision was determined by testing three plasma samples twenty times in one assay.
- Inter-assay precision was determined by testing three plasma samples in twenty assays.

	Intra-Assay Precision			Inter-Assay Precision		
Sample	1	2	3	1	2	3
n	20	20	20	20	20	20
CV (%)	5.0%	3.5%	5.2%	10.3%	9.0%	11.0%
Average CV (%)	4.6%				10.1%	

#### Recovery

Standard Added Value	1 – 10 ng/ml	
Recovery %	89-114%	
Average Recovery %	97%	

#### Linearity

• Plasma and serum samples were serially diluted to test for linearity.

Average Percentage of Expected Value (%)			
Sample Dilution	Plasma	Serum	
1x	94%	90%	
2x	103%	101%	
4x	103%	111%	

### **Cross-Reactivity**

Species	Cross-Reactivity (%)
Canine	>80%
Bovine	None
Equine	10%
Monkey	>80%
Mouse	<10%
Rat	<30%
Swine	>80%
Rabbit	None

• No significant cross-reactivity observed with annexin A1, A2, A3, A4, A8, A10, A11, and A13.

## Troubleshooting

Issue	Causes	Course of Action
	Use of improper components	<ul> <li>Check the expiration date listed before use.</li> <li>Do not interchange components from different lots.</li> </ul>
E	Improper wash step	<ul> <li>Check that the correct wash buffer is being used.</li> <li>Check that all wells are empty after aspiration.</li> <li>Check that the microplate washer is dispensing properly.</li> <li>If washing by pipette, check for proper pipetting technique.</li> </ul>
cisio	Splashing of reagents while loading wells	<ul> <li>Pipette properly in a controlled and careful manner.</li> </ul>
Low Precision	Inconsistent volumes loaded into wells	<ul> <li>Pipette properly in a controlled and careful manner.</li> <li>Check pipette calibration.</li> <li>Check pipette for proper performance.</li> </ul>
	Insufficient mixing of reagent dilutions	<ul> <li>Thoroughly agitate the lyophilized components after reconstitution.</li> <li>Thoroughly mix dilutions.</li> </ul>
	Improperly sealed microplate	<ul> <li>Check the microplate pouch for proper sealing.</li> <li>Check that the microplate pouch has no punctures.</li> <li>Check that three desiccants are inside the microplate pouch prior to sealing.</li> </ul>
High	Microplate was left unattended between steps	• Each step of the procedure should be performed uninterrupted.
ow or ensity	Omission of step Steps performed in incorrect order	<ul> <li>Consult the provided procedure for complete list of steps.</li> <li>Consult the provided procedure for the correct order.</li> </ul>
Unexpectedly Low or High Signal Intensity	Insufficient amount of reagents added to wells	<ul><li>Check pipette calibration.</li><li>Check pipette for proper performance.</li></ul>
Sig	Wash step was skipped	<ul> <li>Consult the provided procedure for all wash steps.</li> </ul>
ue)	Improper wash buffer	<ul> <li>Check that the correct wash buffer is being used.</li> </ul>
Ď	Improper reagent preparation	<ul> <li>Consult reagent preparation section for the correct dilutions of all reagents.</li> </ul>

	Insufficient or prolonged incubation periods	Consult the provided procedure for correct incubation time.
Deficient Standard Curve Fit	Non-optimal sample dilution	<ul> <li>Sandwich ELISA: If samples generate OD values higher than the highest standard point (P1), dilute samples further and repeat the assay.</li> <li>Competitive ELISA: If samples generate OD values lower than the highest standard point (P1), dilute samples further and repeat the assay.</li> <li>User should determine the optimal dilution factor for samples.</li> </ul>
anda	Contamination of reagents	<ul> <li>A new tip must be used for each addition of different samples or reagents during the assay procedure.</li> </ul>
nt Sta	Contents of wells evaporate	<ul> <li>Verify that the sealing film is firmly in place before placing the assay in the incubator or at room temperature.</li> </ul>
Deficie	Improper pipetting	<ul> <li>Pipette properly in a controlled and careful manner.</li> <li>Check pipette calibration.</li> <li>Check pipette for proper performance.</li> </ul>
	Insufficient mixing of reagent dilutions	<ul> <li>Thoroughly agitate the lyophilized components after reconstitution.</li> <li>Thoroughly mix dilutions.</li> </ul>

#### References

- Concha NO, Head JF, Kaetzel MA, Dedman JR, Seaton BA. (Dec 1992). "Annexin V forms calcium-dependent trimeric units on phospholipid vesicles". *FEBS*. 314(2):159-62. doi:10.1016/0014-5793(92)80964-I.
- Sofiadis A, Becker S, Hellman U, Hultin-Rosenberg L, Dinets A, Hulchiy M, Zedenius J, Wallin G, Foukakis T, Höög A, Auer G, Lehtiö J, Larsson C. (April 2012).
   "Proteomic profiling of follicular and papillary thyroid tumors". *Eur J Endocrinol*. 166(4):657-67. doi:10.1530/EJE-11-0856. PMID 22275472.
- (3) Wen, Y, Edelman J L, Kang T, Sachs G. (May 1999). "Lipocortin V may function as a signaling protein for vascular endothelial growth factor receptor-2/Flk-1". *Biochem Biophys Res Commun.* 258(3):713-21. doi:10.1006/bbrc.1999.0678. ISSN 0006-291X. PMID 10329451.
- (4) Cardó-Vila M, Arap W, Pasqualini R. (May 2003). "Alpha v beta 5 integrindependent programmed cell death triggered by a peptide mimic of annexin V". *Mol Cell*. 11(5):1151-62. doi:10.1016/S1097-2765(03)00138-2. ISSN 1097-2765. PMID 12769841.

Version 1.4R